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
REF 053 900

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D053900I2022-04I007 ML I Rev. 01

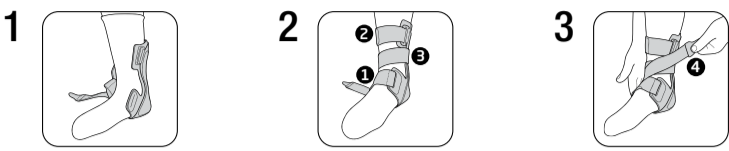DE BORT TaloXpress

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EN BORT TaloXpress

Many thanks for placing your trust in a medical device from BORT GmbH. Please read the existing instructions for use carefully. If you have any questions, please refer to your physician or the specialist retailer from whom you purchased this medical device.

Intended purpose

This medical device is an adjustable ankle brace for the stabilisation of the upper and lower ankle joint on one plane. The brace consists of elastic material and stabilisation elements made of rigid material.

Indications

Conservative treatment of ankle distortions and ligament ruptures, post-operative protection after ligament suture/reconstruction, immobilisation in the case of inflammatory conditions, rheumatism, arthritis, tendinitis, supination/pronation trauma

Contraindications

Thrombosis risk, extreme varicosis, peripheral arterial occlusive disease (PAOD), lymph drainage disorders, also unclear soft tissue swellings distal to the aid positioned, sensory and circulatory disorders of the region of the body treated, skin diseases in the part of the body treated

Application risks/important notes

- This medical device is a prescribable product. Discuss the use and duration with your treating physician. The expert staff from whom you have received the medical device will select the appropriate size and instruct you regarding its use.
- remove the medical device prior to radiological examinations
- If wearing at night is prescribed, avoid negative influence on the circulatory system
- in case of numbness, loosen the medical device and remove it if necessary
- in case of persistent complaints, consult the physician or specialist retailer
- use the medical device in accordance with therapeutic needs
- only use other products simultaneously after consultation with your physician
- do not make any changes to the product
- do not wear it on open wounds
- do not use in case of intolerance of one of the materials used
- no re-use – this medical aid is intended for treating one patient
- whilst wearing the brace, please neither use any creams nor ointments on or around it as they can destroy the material

Fitting instructions

Important: The BORT TaloXpress is worn over the socks, not directly on the skin.

Putting the brace on:

Open all the brace's fasteners. For this purpose, unloop the belt straps. If required, the strap length can be set by shortening the straps.

1 Place your foot into the brace, here the plastic shell must be positioned on the inner side of the leg. The lower leg should be positioned vertical to the ground.

2 Close the fastener on the metatarsal first, it is marked with a "1". Then, fix the fastener on the upper brace edge which is marked with the number "2". In the next step, close fastener "3" on the blue Velcro surface.

3 Finally, fastener "4" is used to set the stabilisation degree. For this purpose, guide the sealing tape diagonally over the instep and fix it to the blue Velcro surface under slight tension.

To improve relief, the outer edge of the foot can be lifted slightly whilst the diagonal restraint is positioned.

If required, the locking belts can be readjusted after the brace has been positioned.

Removal

To remove the BORT TaloXpress, open the belts in the reverse order.

Tip: Close the Velcro fasteners in order to avoid early wear to the Velcro surface.

Material composition

Polyamide (PA), polyester (PES), polyurethane (PUR)

The sewn in textile label provides the precise material composition.

Cleaning information

- Hand-washing Do not bleach Do not dry in a tumble dryer
- Do not iron Do not clean chemically
- Do not use fabric conditioner. Stretch back into shape and dry in the fresh air.
- Close the Velcro fastenings to avoid damage to other laundry items.

Guarantee

The legal regulations of the country in which you acquired the product apply to the purchased product. Please contact your specialist retailer if you suspect a warranty claim. Please clean the product before submitting a warranty claim. If the enclosed instructions for use have not been properly observed, the warranty may be impaired or cancelled. The warranty does not cover use of the product inappropriate for the indication, non-observance of application risks, instructions and unauthorised modifications to the product.

Useful life/Lifetime of the product

The lifetime of the medical device is determined by its natural wear and tear if treated appropriately and as recommended.

Duty of notification

If a serious worsening of the state of health occurs when using the medical device, you are obliged to notify the specialist dealer or us as the manufacturer and the Medicines & Healthcare products Regulatory Agency (MHRA). You can find our contact information in these instructions for use. You can find the contact information for the national authority responsible for your country under the following link: www.bort.com/md-eu-contact.

Disposal

Upon the termination of use, the product must be disposed of in accordance with the corresponding local requirements.

Declaration of conformity

We confirm that this device conforms with the requirements of REGULATION (EU) 2017/745 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL. You can find the current declaration of conformity under the following link: www.bort.com/conformity

Status: 04.2022

Medical device | Single patient – multiple use

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